

ri-medic



ri-medic clinical-grade blood pressure monitor.

One touch automated and manual blood pressure. World class, clinical grade technology for fast and reliable readings with full range of blood pressure cuff sizes. Available with mobile stand, desk top or wall mount.

- World class, clinical grade NIBP.
- BP independently validated to British Hypertension Society, ESH and AAMI SP10 specifications for adults and children (from 3 years onwards).
- Manual SphygMode with auto-deflate provides consistent patient BP measurement and allows clinicians to verify automated readings without removing the cuff.
- Designed for durability, comes with a 2-year warranty.
- Easy to use one touch operation.
- Fast readings, typically 30 – 40 seconds.
- Mobile versions (tabletop and on mobile stand) with rechargeable internal battery.
- Storage basket on the reverse side for mobile versions and for wall-mount (included), optional basket attached to pole.
- Optional: Bluetooth® wireless EMR/EHR interface with Management Software (510 K listed).

ri-medice
clinical-grade blood pressure monitor.

Technical data

Automated BP Specification:

BP Measurement: Oscillometric
Measurement Range: Systolic: 60-270 mmHg
 Diastolic: 30-170 mmHg
 Heart Rate: 30-200 bpm
 MAP: 40-203 mmHg
BP Read Time: Typically: 30-40 seconds
Accuracy: BHS (A/A), ESH, AAMI SP10:2002
Heart Rate Accuracy: +- 2% or +- 3bpm, whichever is greater

Manual SphygMode BP Specification:

Deflate Rate: 3 mmHg/sec (Meets AHA Recommendations)

Cuff Specifications:

Tube length 70cm/27.56 in (cuff only), total tube length 310cm/122.05 in

Rechargeable internal battery for desk top and mobile version:

200 measurements within a 12-hour period, 8-12 hours per charge.

Mains adapter:

Input: 110-240V AC, 50-60Hz, 0.6A
 Country-specific interchangeable socket plugs available.

ri-medice table top:

Dimensions: W 18cm/7.09 in x L 23cm/9.06 in x H 25cm/9.84 in
Weight: 2200 g / 4.85 lb

ri-medice mobile stand:

- Height adjustable from 75 cm / 29.53 in to 120 cm / 47.24 in
- Diameter 53 cm / 20.87 in
- Sturdy mobile base made of fibreglass and chrome plated pole.
- Easy to move due to low friction castors.
- Optional basket attached to pole.
 Small: L 25cm/9.84in x W 9cm/3.54in x H 20cm/7.87in
 Large: L 23cm/9.06in x W 23cm/9.06in x H 25cm/9.84in
- Power supply attached to the pole.

Ordering Information:

Tabletop with adult cuff and adult obese cuff, latex free	Art.No.: 1780
Tabletop with adult cuff and adult obese cuff, latex free, Bluetooth®	Art.No.: 1781
Mobile version with adult cuff and adult obese cuff, latex free	Art.No.: 1782
Mobile version with adult cuff and adult obese cuff, latex free, Bluetooth®	Art.No.: 1783
Wall mount with adult cuff and adult obese cuff, latex free	Art.No.: 1784
Wall mount with adult cuff and adult obese cuff, latex free, Bluetooth®	Art.No.: 1785
Mobile stand only (for tabletop)	Art.No.: 10395
Small basket for mobile stand	Art.No.: 10392
Large basket for mobile stand	Art.No.: 10393
Velcro cuff adult, latex free	LF107
Velcro cuff adult obese, latex free	LF122
Velcro cuff adult small, latex free	LF120
Velcro cuff children, latex free	LF129




ri-medice table top



ri-medice wall mount



ri-medice mobile stand

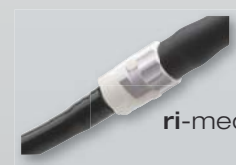

ri-medice cuff quick connect